



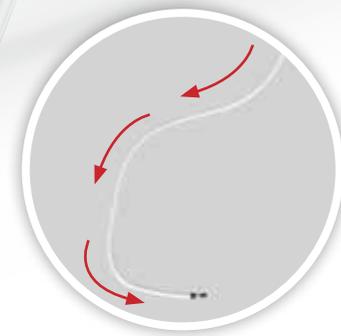

CARDIAC INTERVENTION

# Worley™ Advanced RCSG

## Right Sided Coronary Sinus Guide

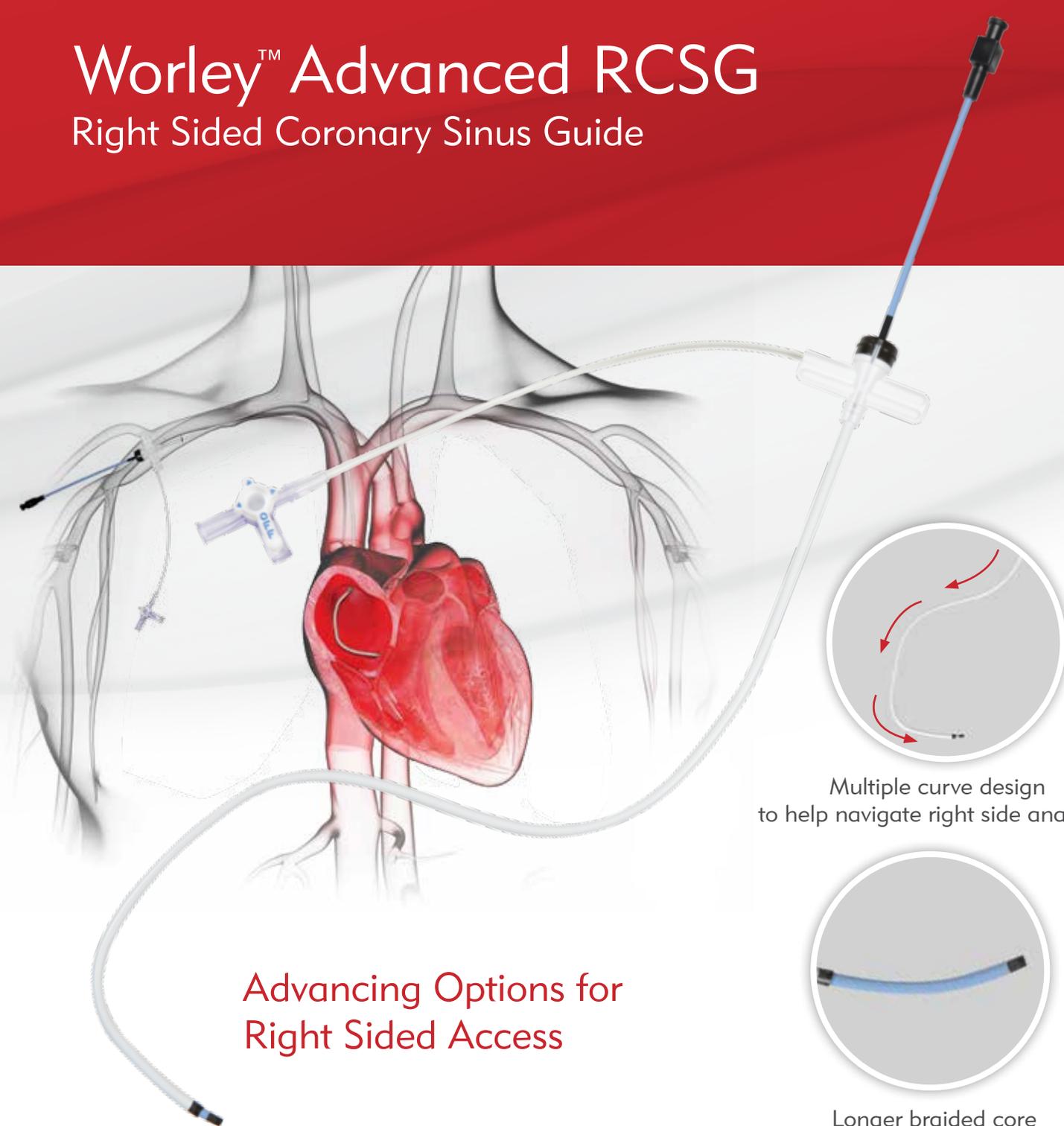

Multiple curve design to help navigate right side anatomy



Longer braided core for extended reach

Advancing Options for Right Sided Access

# Worley™ Advanced RCSG

Right Sided Coronary Sinus Guide

## PROMOTING RELIABLE ACCESS & ACCURATE DELIVERY

To meet the challenges of right sided access during a cardiac resynchronization therapy (CRT), you need an advanced LV delivery system that gives you options to match the tortuosity of your patient's anatomy. The Worley Advanced Right Sided CSG is designed to provide you the options you need to reliably and accurately place your LV leads at the target location.

## PRODUCT FEATURES

**MULTIPLE SHEATH CURVES** to provide back wall support while navigating the right subclavian venous anatomy

**BRAIDED CORE** with its longer variable shape design provides control and flexibility

**DEPTH POSITION MARKER** for alignment of the distal tips of the braided core and the RCSG

**INTEGRATED PEEL-AWAY HEMOSTASIS VALVE** handles like a standard peel-away introducer to reduce risk of blood loss and air embolism


**SIDE PORT** allows open access infusion and contrast injection during implantation

**CRUSH RESISTANT PEBAX®\* SHEATH** with braided core reduces kinking and facilitates easy removal

**9 French ID** provides ample room to perform various advanced CRT implantation techniques

**SOFT RADIOPAQUE TIP** reduces risk of vessel wall trauma and improves fluoroscopic visualisation

## VARIABLE SHAPE DESIGN

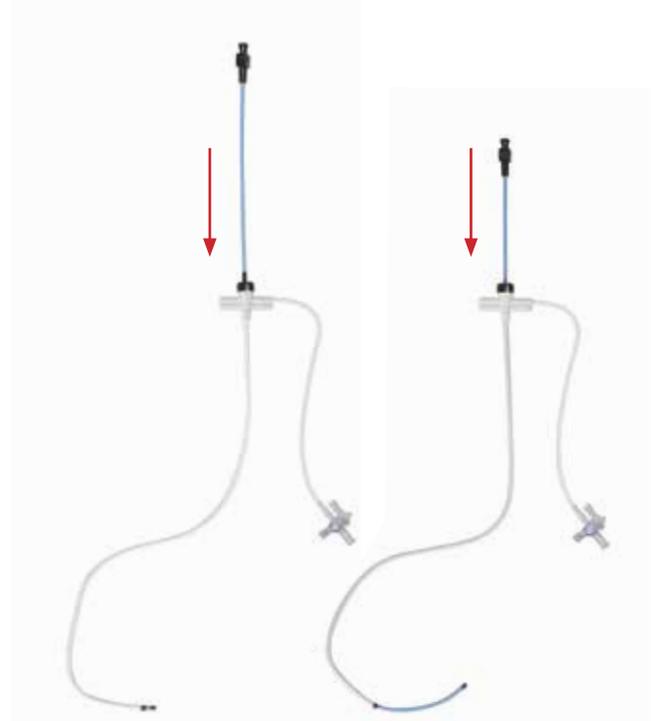

Part Number	French Size	Length
WOR-CSG-BR1-09	9F	40 cm

### Kit Contents:

- Splittable Sheath
- Dilator
- Braided Core
- Guidewires:
  - 0.035" x 135 cm
  - 0.032" x 180 cm
- Transvalvular Insertion Tool
- 18G 7 cm Needle
- 12 mL Syringe



Understand. Innovate. Deliver.™

Merit Medical Systems, Inc.  
1600 West Merit Parkway  
South Jordan, Utah 84095  
1.801.208.4300  
1.800.35.MERIT

Merit.com

Merit Medical Europe, Middle East, & Africa (EMEA)  
Amerikalaan 42, 6199 AE Maastricht-Airport  
The Netherlands  
+31 43 358 82 22

Merit Medical Ireland Ltd.  
Parkmore Business Park West  
Galway, Ireland  
+353 (0) 91 703 733

Austria  
0800 295 374

Belgium  
0800 72 906 (Dutch)  
0800 73 172 (Français)

Denmark  
80 88 00 24

Finland  
0800 770 586

France  
0800 91 60 30

Germany  
0800 182 0871

Ireland (Republic)  
1800 553 163

Italy  
800 897 005

Luxembourg  
8002 25 22

Netherlands  
0800 022 81 84

Norway  
800 11629

Portugal  
308 801 034

Russia  
+7 495 221 89 02

Spain  
+34 911238406

Sweden  
020 792 445

Switzerland  
(Deutsch)  
+41 225180252  
(Français)  
+41 225948000  
(Italiano)  
+41 225180035  
UK  
0800 973 115